

# Model FDL 115 | Safety drying chambers for limited quantities of solvent

The BINDER FDL 115 with silicone and dust-free inner chamber, as well as symmetric airflow dries specimens containing solvents and thus this safety drying oven meets all EN 1539 requirements.

# BENEFITS

- Safety concept meets the DIN EN 1539-2010-08 standard
- Short drying time with high temperature accuracy and defined ventilation
- Silicone free and dust-free inner chamber
- Replaceable fresh air cartridges are easy to change and ensure constant fresh air quality





Model 115

## **MAIN FEATURES**

- Temperature range: +10 °C above ambient temperature to +300 °C
- APT.line™ preheating chamber technology
- Controller with time-segment and real-time programming
- 2 chrome-plated racks
- Fresh-air monitoring with audible and visual alarm and automatic deactivation of heating
- Replaceable filter cartridge, Class M6

- Class 2 independent adjustable temperature safety device (DIN 12880) with visual alarm
- Computer interface: RS 422
- Standard door gasket made from solvent-resistant FKM for use at up to 200  $^{\circ}\mathrm{C}$
- Optional door gasket made from silicone with limited solvent resistance for use at up to 300 °C

## **ORDERING INFORMATION**

Interior volume [L]	Voltage	Option model	Version	ArtNo.
115	230 V 1~ ph 50/60 Hz	Standard	FDL115-230V	9010-0292




# **TECHNICAL DATA**

Description	FDL115-230V <sup>1</sup>
Article Number	9010-0292
Performance Data Temperature	
Temperature range +10 °C above ambient temperature to [°C]	300
Temperature variation at 150 °C [± K]	2.5
Temperature fluctuation [± K]	0.8
Heating-up time to 150 °C [min]	20
Recovery time after 30 seconds door open at 150 °C [min]	12
Air change data	
Air circulation (approx.) [x/min]	20
Volumetric flow rate of exhaust air acc. to EN 1539 at 50 °C [ca. L/min]	400
Solvent quantity	
Highest permitted solvent quantity (at T-180 °C, M-100g/mol, U-40g/m3, K=0,5) [g]	6.65
Electrical data	
Rated Voltage [V]	230
Power frequency [Hz]	50/60
Nominal power [kW]	2.9
Phase (Nominal voltage)	1~
Measures	
Interior volume [L]	115
Net weight of the unit (empty) [kg]	90
Permitted load [kg]	60
Load per rack [kg]	30
Wall clearance back [mm]	160
Wall clearance sidewise [mm]	100
Internal Dimensions	
Width [mm]	600
Height [mm]	435
Depth [mm]	435
Housing dimensions not incl. fittings and connections	
Width net [mm]	830
Height net [mm]	805
Depth net [mm]	685
Environment-specific data	
Sound-pressure level [dB(A)]	57
Fixtures	
Number of shelves (std./max.)	2/5

1 All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10 %. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. Technical data refers to 100 % fan speed. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.

## DIMENSIONS incl. fittings and connections [mm]

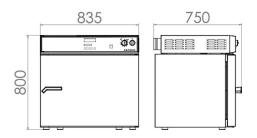



## **OPTIONS**

Designation	Description	*	ArtNo.
	back		
	10 mm	01	8012-1284
	left		
Access port with silicone	10 mm	01	8012-1277
plug	30 mm	01	8012-1360
	right		
	10 mm	01	8012-1270
	30 mm	01	8012-1355
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test temperature	-	8012-1117
	for temperature, measurement in center of chamber at specified temperature	-	8012-1136
Calibration certificate,	temperature measurement incl. certificate and 27 measuring points at specified temperature	-	8012-1595
temperature	temperature measurement incl. certificate, 15-18 measuring points at specified temperature	-	8012-1574
	temperature measurement incl. certificate, 9 measuring points at specified temperature	-	8012-1553
Coil coating door flap	door flap for coil coating tests	-	8012-0179
Door gasket	made from silicone, resistant to high temperatures for operating temperatures above 200 °C. Note: The gasket wears out more quickly in temperatures above 250 °C	-	8012-0499
Door lock	lockable door handle	-	8012-1248
Measurement of air exchange rate	in accordance with ASTM D5374, definition and protocol according to ambient temperature	_	8012-1196
* Notos - Cas last same			

\* Notes > See last page

# ACCESSORIES

Designation	Description	*	ArtNo.
APT-COM™ 4 BASIC-	for simple logging and documentation requirements with up to 5 networked units.		
Edition	version 4, BASIC edition	-	9053-0039
APT-COM™ 4 GLP-	for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11.		
dition	version 4, GLP edition	-	9053-0042
PT-COM™ 4	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.		
PROFESSIONAL-Edition	version 4, PROFESSIONAL edition	-	9053-0040
H-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	-	1002-0016
	IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature, CO <sub>2</sub> , O <sub>2</sub> , pressure, depending on unit		
	Digital in PDF format	-	7057-0001
	Hard copy inside folder	-	7007-0001
Qualification documents	IQ/OQ/PQ documents – supporting documents for validation performed by customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature, $CO_2$ , $O_2$ – or pressure, depending on unit		
	Digital in PDF format	-	7057-000
	Hard copy inside folder	-	7007-000
Rack	chrome plated	-	8012-2034
	stainless steel	-	8012-2048
helf, perforated	Stainless steel	-	6004-003
able on castors	stable cart, casters with locking brakes, dimensions: W 1,000 x D 800 x H 780 mm	_	9051-0018

\* Notes > See last page

# SERVICES



Designation	Description	*	ArtNo.
Installation services			
Unit commissioning	Connect the unit to the customer-side connections (electricity, water, wastewater, gas), basic functions check, brief operating instructions. (excl.: unpacking, setup, controller instructions, programming, installation work)	05	DL10-0100
Unit instructions	Instruction regarding operating principle and basic functions of the unit, operation of the control electronics including programming	05	DL10-050
Maintenance contracts			
Naintenance services			
3RONZE 3-year maintenance contract	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts	05	DL20-071
GOLD 3-year maintenance contract	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts, testing of all key functions, replacement of wear parts, calibration of one test temperature specified by the user in the center of the usable space, including certificate	05	DL20-091
Maintenance	One-off maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Calibration of a test temperature specified by the user in center of usable space without certificate	05	DL20-020
SILVER 3-year maintenance contract	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts, testing of all key functions, calibration of one test temperature specified by the user in the center of the usable space, without certificate	05	DL20-081
Calibration services			
Calibration certificate, temperature	Calibration of one (1) test temperature specified by the user in center of chamber, including certificate	03, 04, 05	DL30-010
	Extension of calibration of one (1) additional test temperature specified by the user in the center of the usable space, including certificate	03, 04, 05	DL30-010
Measurement of air exchange rate	including certificate (in accordance with ASTM D5374)	03, 04, 05	DL33-000
emperature measurement, 18 neasuring points	Temperature measurement with 18 measuring points with a set value specified by the user, including certificate	03, 04, 05	DL30-0118
emperature measurement, 27 neasuring points	Temperature measurement with 27 measuring points with a set value specified by the user, including certificate	03, 04, 05	DL30-012
emperature measurement, 9 neasuring points	Temperature measurement with 9 measuring points with a set value specified by the user, including certificate	03, 04, 05	DL30-010
alidation services			
xecution of IQ/OQ	Execution of IQ/OQ in accordance with qualification folder	05	DL40-010
execution of IQ/OQ/PQ	Execution of IQ/OQ/PQ in accordance with qualification folder	05	DL44-050
Narranty service			
year warranty extension	The warranty is extended by 1 year from the delivery date, wear parts are excluded	-	DL50-001
Notes > See last nage			

\* Notes > See last page



#### NOTES

- Condensation may occur in the area around the access port. Access ports may be placed in custom locations for an additional charge. UL mark is not granted when this option is used. 01
- 02
- Sensor calibration is performed in an accredited calibration laboratory. 03
- Calibration is performed according to the BINDER factory standard. 04
- 05 Quoted prices do not include travel costs. Please refer to the chapter on BINDER Service for travel costs for your region. Quoted prices for services performed in Switzerland do not include a country-specific added fee (available on request).

#### **BINDER GmbH**

Tuttlingen, Germany TEL +49 7462 2005 0 FAX +49 7462 2005 100 info@binder-world.com www.binder-world.com

#### BINDER Asia Pacific (Hong Kong) Ltd.

Kowloon, Hong Kong, P.R. China TEL +852 39070500 FAX +852 39070507 asia@binder-world.com www.binder-world.com

#### **BINDER Environmental Testing**

Equipment (Shanghai) Co., Ltd. Shanghai, P.R. China TEL +86 21 685 808 25 FAX +86 21 685 808 29 china@binder-world.com www.binder-world.com

### **Representative Office for CIS**

Moscow, Russia TEL +7 495 988 15 16 FAX +7 495 988 15 17 russia@binder-world.com www.binder-world.com

#### **BINDER Inc.**

Bohemia, NY, USA TEL +1 631 224 4340 FAX +1 631 224 4354 usa@binder-world.com www.binder-world.us